


UFII Master Kit

Kit Code | UF(M) 06

- Kit composed of drill with outstanding cutting force and durability, etc.
- UFII Narrow Ø3.0 Ø3.3/ UFII Regular Ø3.8 Ø4.0 Ø4.5 Ø5.0 Fixture

Unit mm | Drill Scale 1 : 0.9

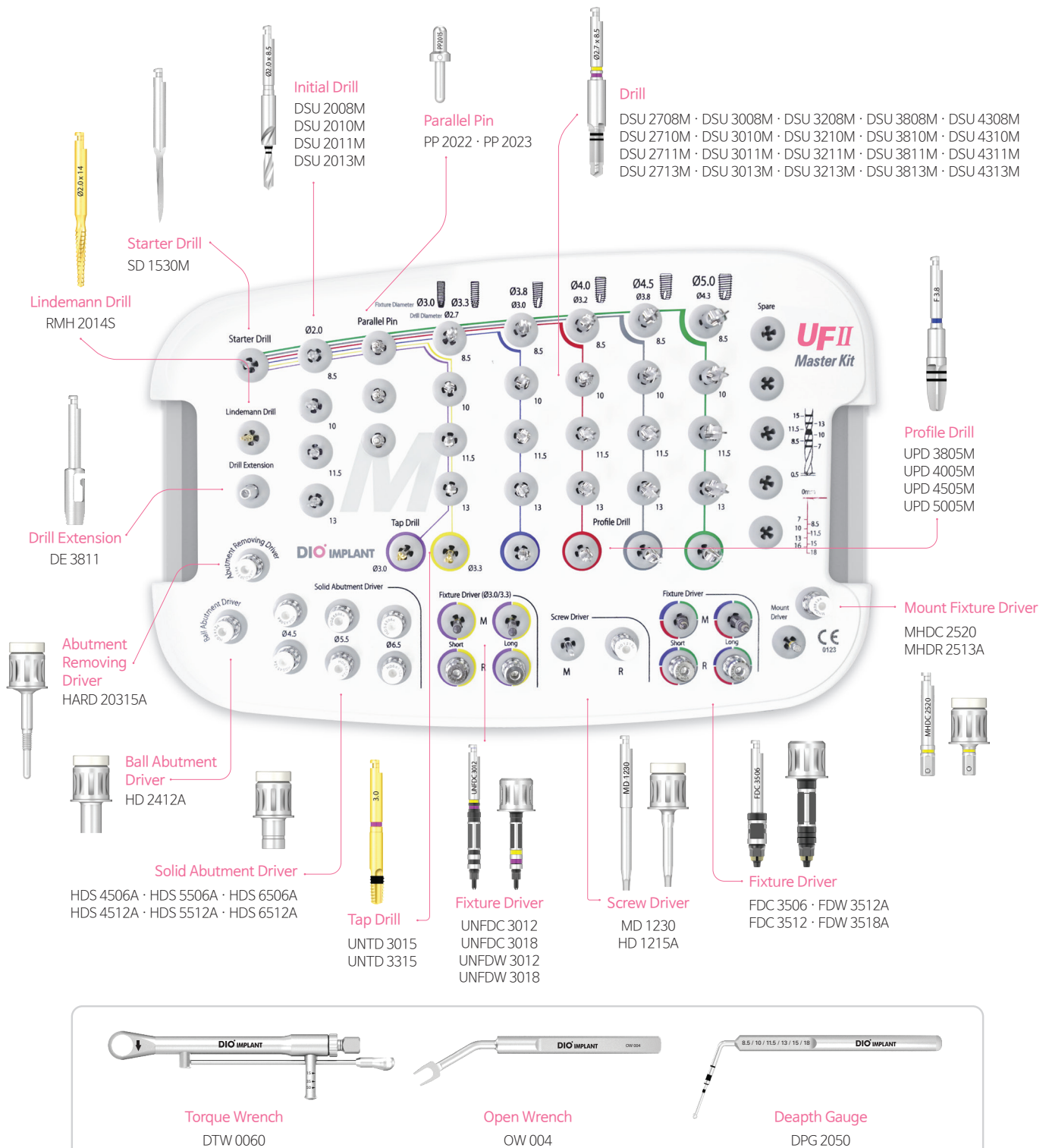

Ver A



UFII Surgical Kit Kit Code | UF(M) 07

- Kit composed of drill with outstanding cutting force and durability, etc.
- UFII Narrow Ø3.0 Ø3.3/ UFII Regular Ø3.8 Ø4.0 Ø4.5 Ø5.0 Ø5.5 Fixture

Unit mm | Drill Scale 1 : 0.9

IMPLANT SYSTEM

UFII IMPLANT SYSTEM | UFII Master Surgical Kit

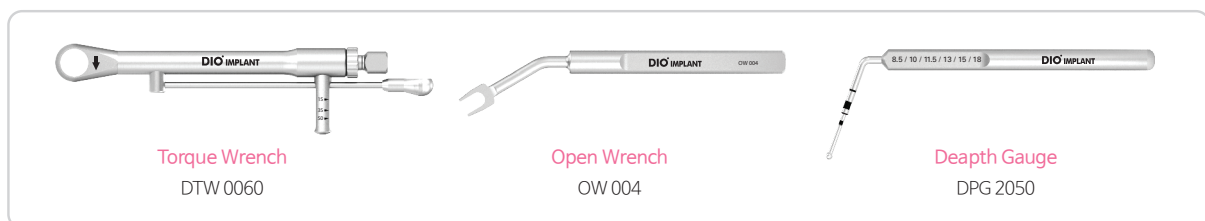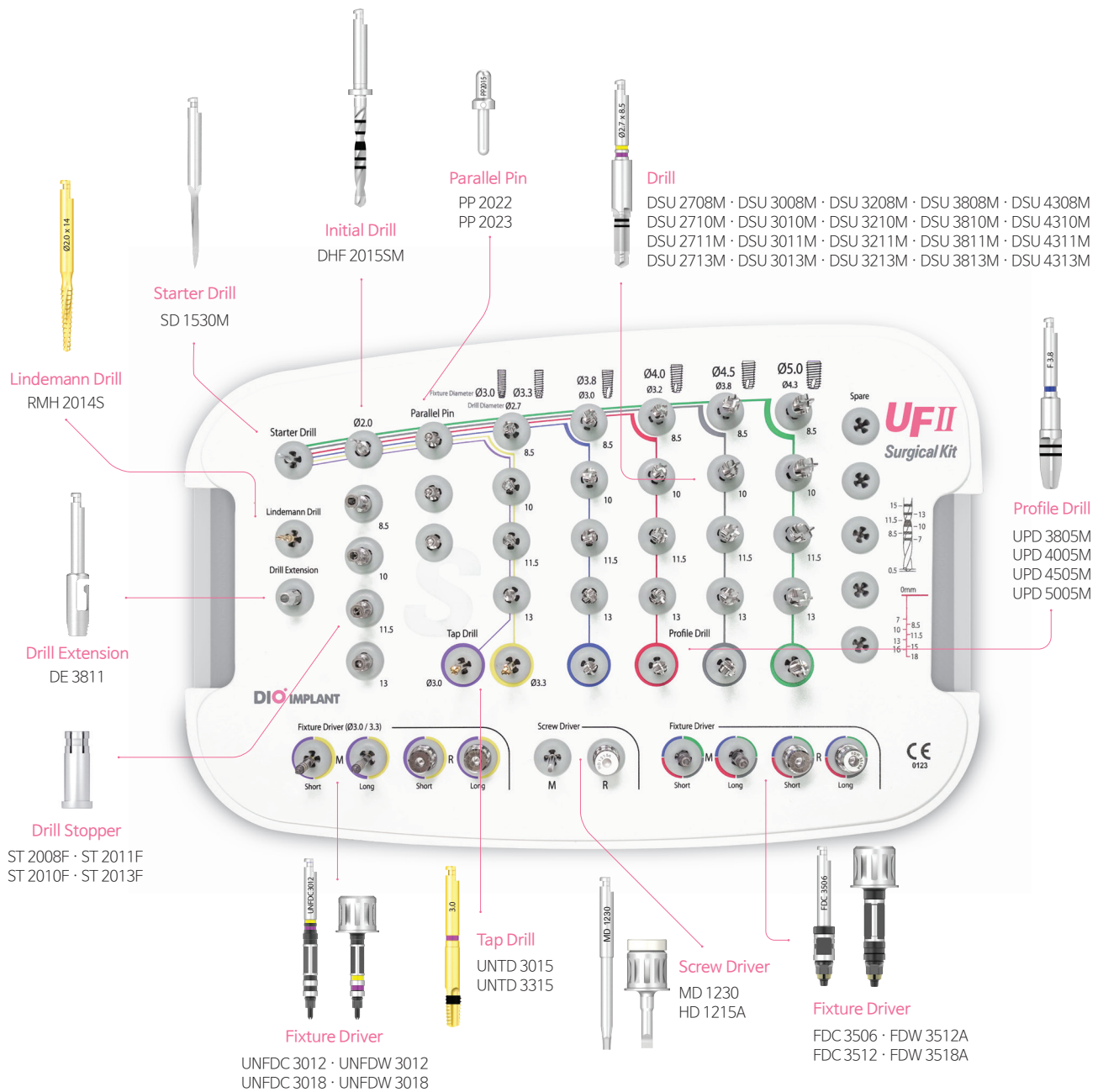


Ver B

Surgical Tool


Starter Drill

K	
Code	SD 1530M

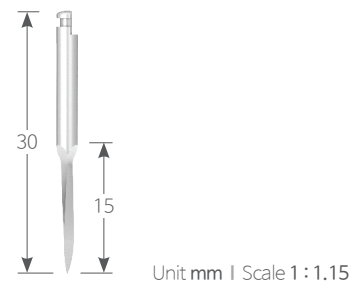


Initial Drill

K D Ø2.0 · Only use UFII Master Kit		Unit mm Scale 1 : 1		
Code	DSU 2008M	DSU 2010M	DSU 2011M	DSU 2013M

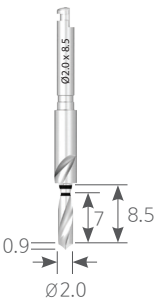

K D Ø2.0 · Only use UFII Surgical Kit		Unit mm Scale 1 : 1.2	
Code	DHF 2015SM		

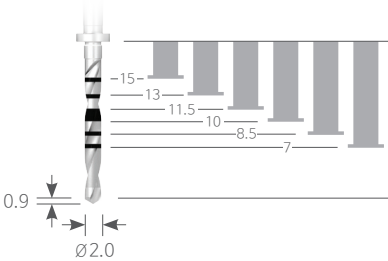

Drill stopper

K D Ø2.0 · Only use UFII Surgical Kit		Unit mm Scale 1 : 1		
Code	ST 2008F	ST 2010F	ST 2011F	ST 2013F



Parallel Pin


● Sold separately

K	Unit mm Scale 1 : 1				
	Code	PP 1223	PP 1523	PP 2015	PP 2022 PP 2023

Path Pin

● Sold separately


K	Code	DPP 3512
---	------	----------

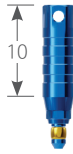

Lindemann Drill

Unit mm | Scale 1 : 1.2

K	D Ø2.0	RMH 2014S
---	--------	-----------

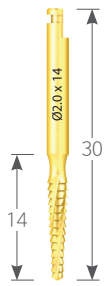

Positioning Guide

● Sold separately

K	D Ø2.0	PG 0060
---	--------	---------



Angled Path Pin

● Sold separately

Unit mm | Scale 1 : 1.15

K	Code	DAP 4515A DAP 4515B
---	------	------------------------



Type

Hex A

Hex B

Surgical Tool

Drill



● Sold separately

K D Ø2.7

Unit mm | Scale 1 : 1

Code	DSU 2708M	DSU 2710M	DSU 2711M	DSU 2713M	DSU 2715M	DSU 2718M
						
Length	8.5	10	11.5	13	15	18

K D Ø3.0

Unit mm | Scale 1 : 1

Code	DSU 3008M	DSU 3010M	DSU 3011M	DSU 3013M	DSU 3015M	DSU 3018M
						
Length	8.5	10	11.5	13	15	18

K D Ø3.2

Unit mm | Scale 1 : 1

Code	DSU3208M	DSU 3210M	DSU 3211M	DSU 3213M	DSU 3215M	DSU 3218M
						
Length	8.5	10	11.5	13	15	18

K D Ø3.8

Unit mm | Scale 1 : 1

Code	DSU 3808M	DSU 3810M	DSU 3811M	DSU 3813M	DSU 3815M	DSU 3818M
						
Length	8.5	10	11.5	13	15	18

K D Ø4.3

Unit mm | Scale 1 : 1

Code	DSU4308M	DSU 4310M	DSU 4311M	DSU 4313M	DSU 4315M	DSU 4318M
						
Length	8.5	10	11.5	13	15	18

K D Ø4.8

Unit mm | Scale 1 : 1

Code	DSU 4808M	DSU 4810M	DSU 4811M	DSU 4813M	DSU 4815M	DSU 4818M
						
Length	8.5	10	11.5	13	15	18

Surgical Tool


Profile Drill

K

Sold separately

Unit mm | Scale 1 : 1

Fixture Size	Ø3.8	Ø4.0	Ø4.5	Ø5.0	Ø5.5
Code	UPD 3805M	UPD 4005M	UPD 4505M	UPD 5005M	UPD 5505M

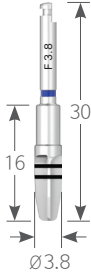

Ø3.8



Ø4.0



Ø4.5



Ø5.0



Ø5.5

Color

●

●

●

●

○

Tap Drill

K

Unit mm | Scale 1 : 1

Fixture Size	Ø3.0	Ø3.3	Ø3.8	Ø4.0	Ø4.5	Ø5.0	Ø5.5
Code	UNTD 3015	UNTD 3315	UT (II) 3815	UT (II) 4015	UT (II) 4515	UT (II) 5015	UT (II) 5515



Ø3.0



Ø3.3



Ø3.8



Ø4.0



Ø4.5



Ø5.0



Ø5.5

Color

●

●

●

●

●

●

○

Drill Extension

K

Unit mm | Scale 1 : 1.2




DE 3811

Code	DE 3811
------	---------

Fixture Driver Narrow

● Sold separately

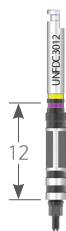Machine		Unit mm Scale 1 : 1		
Code	UNFDC 3006	UNFDC 3012	UNFDC 3018	
				
Ratchet				
Code	UNFDW 3006	UNFDW 3012	UNFDW 3018	

Fixture Driver Regular / Wide

● Sold separately

K



Machine

Unit mm | Scale 1 : 1

Code	FDC 3506	FDC 3512	FDC 3518
			

K

Ratchet

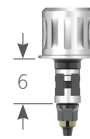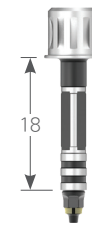Code	FDW 3506A	FDW 3512A	FDW 3518A
			

Surgical Tool

Screw Driver Machine

● Sold separately

K Slot 0.5


Unit mm | Scale 1 : 1

Code	MD 0522	MD 0530
------	---------	---------



K Hex 0.9

Code	MD 0922	MD 0930
------	---------	---------

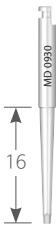


K Hex 1.2

Code	MD 1219	MD 1222	MD 1230	MD 1234	MD 1239
------	---------	---------	---------	---------	---------



K Torx 1.7

Code	MD 1719	MD 1722	MD 1730
------	---------	---------	---------


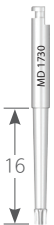


Screw Driver Ratchet




● Sold separately

K Slot 0.5







Unit mm | Scale 1 : 1

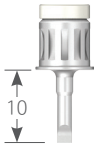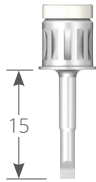Code	HD 0510A	HD 0515A
		






K Hex 0.9

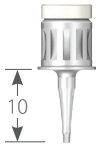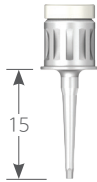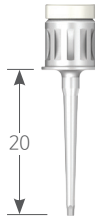Code	HD 0910A	HD 0915A	HD 0920A
			

K Hex 1.2

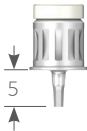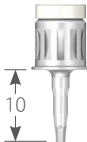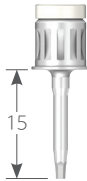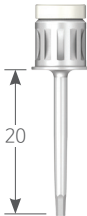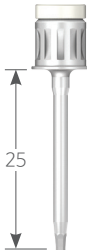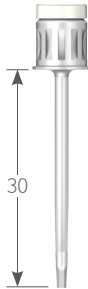Code	HD 1205A	HD 1210A	HD 1215A	HD 1220A	HD 1225A	HD 1230A
						

K Torx 1.7

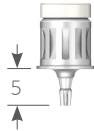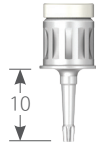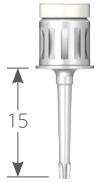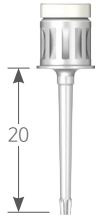
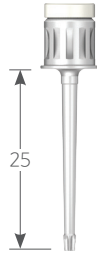
Code	HD 1705A	HD 1710A	HD 1715A	HD 1720A	HD 1725A
					

Surgical Tool

Fixture Mount Driver

· Only use UF II Master Kit



K Hex 2.5		Unit mm Scale 1 : 1	
Type	Machine		Ratchet
Code	MHDC 2520	MHDC 2525	MHDR 2513A MHDR 2518A

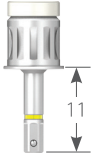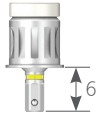

Solid Abutment Driver

· Only use UF II Master Kit

K		Unit mm Scale 1 : 1.2	
Length	6		12

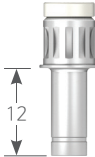
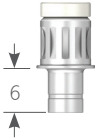


Diameter Ø4.5	HDS 4506A	HDS 4512A
Ø5.5	HDS 5506A	HDS 5512A
Ø6.5	HDS 6506A	HDS 6512A

Ball Abutment Driver



· Only use UF II Master Kit

● Sold separately

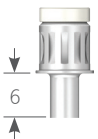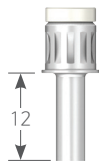K	Unit mm Scale 1 : 1	
	Code	HD 2406A HD 2412A
		


Abutment Removing Driver

· Only use UF II Master Kit

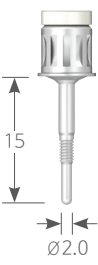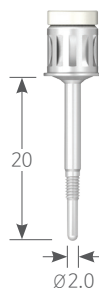K	Unit mm Scale 1 : 1	
	Code	HARD 20315A HARD 20320A
		

Surgical Tool

Depth Gauge

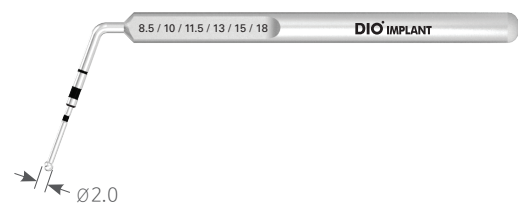

K D Ø2.0

Unit mm | Scale 1 : 1

Code

DPG 2050

Open Wrench

- Used to remove Mount if the bone is soft.
- 30 ° Neck angle for convenience in intraoral insertion.



K

Unit mm | Scale 1 : 1

Code

OW 004

Torque Wrench



K


Unit mm | Scale 1 : 1

Code


DTW 0060

Surgical Protocol


Crestal Placement

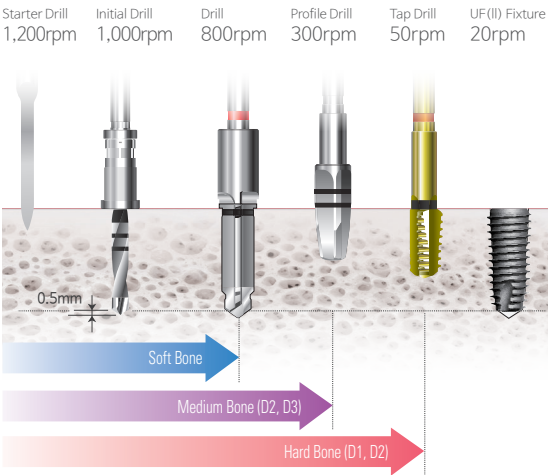

1mm Subcrestal Placement

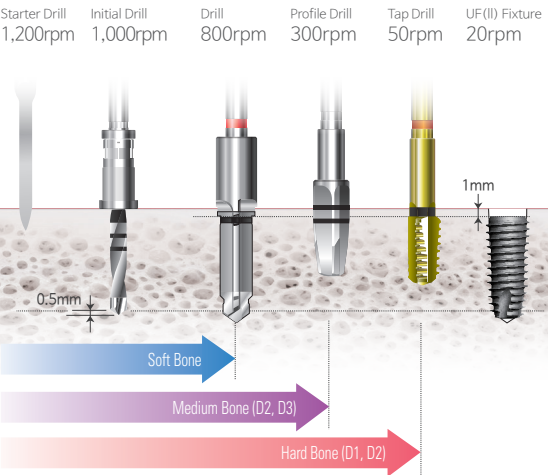


Drill
Using 8.5mm Drill with laser marking to place 7mm implant

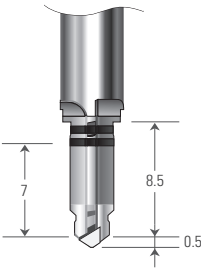


Profile Drill

Drilling depth can effect initial stability.

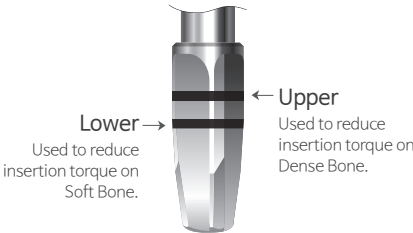

Tap Drill


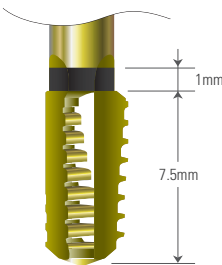


UF II Narrow Kit

Kit Code | UF(M) 08


- Kit composed of drill with outstanding cutting force and durability, etc.
- UF II Narrow Ø3.0 Ø3.3 Fixture

Unit mm | Drill Scale 1 : 0.9

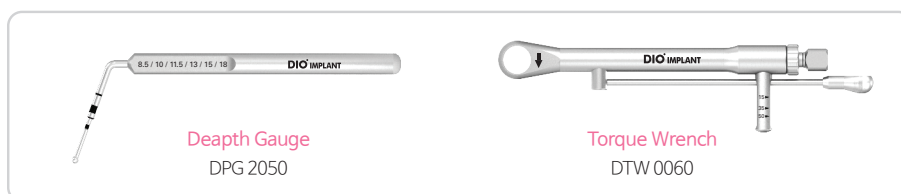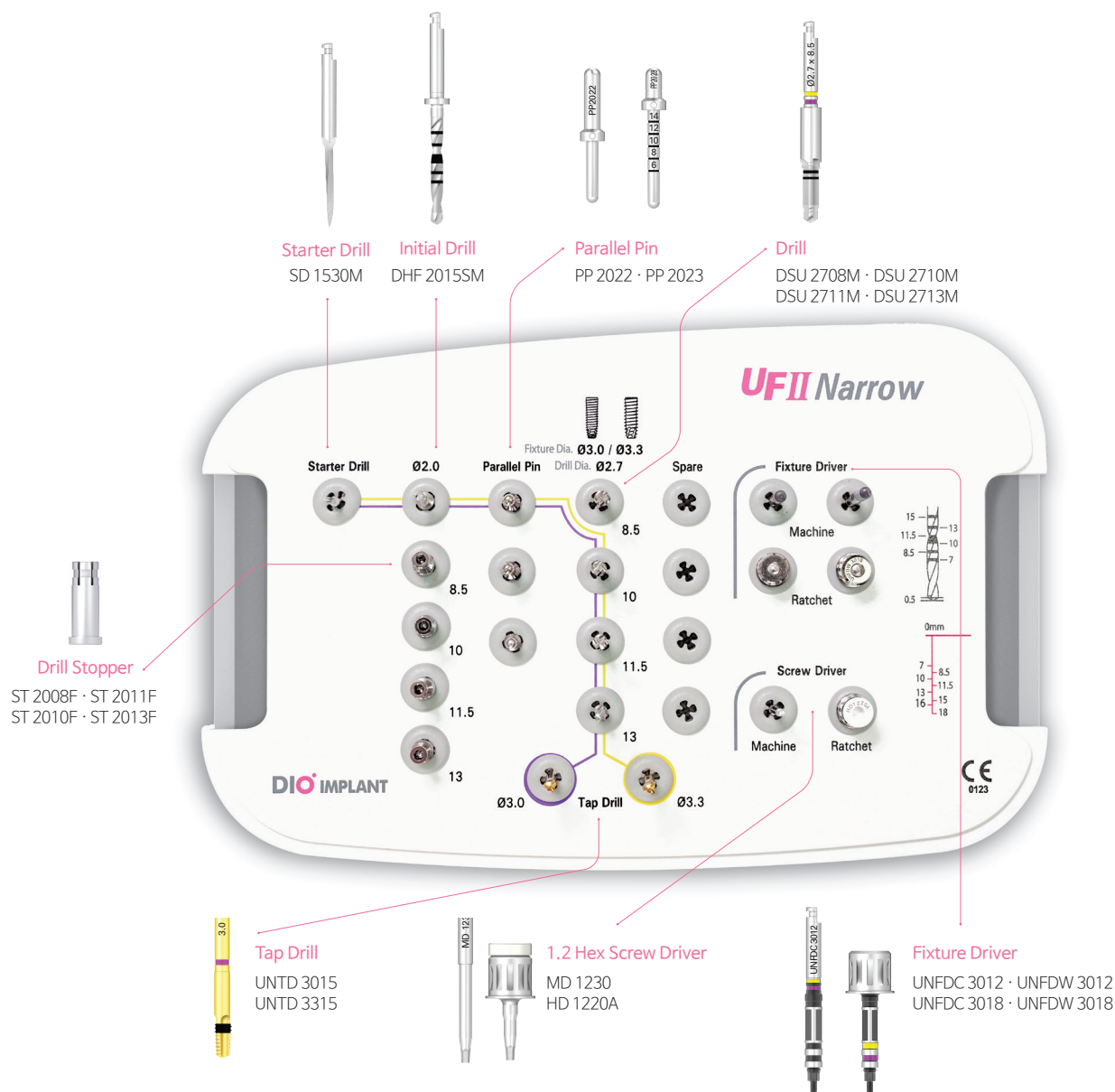


Surgical Tool

Starter Drill

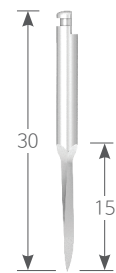

K	Unit mm Scale 1 : 1.2	
Code	SD 1530M	

Initial Drill

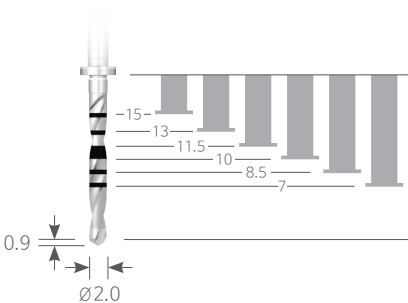



K D Ø2.0	Unit mm Scale 1 : 1.2	
Code	DHF 2015SM	

Drill stopper

K D Ø2.0	Unit mm Scale 1 : 1				
Code	ST 2015F	ST 2013F	ST 2011F	ST 2010F	ST 2008F
					


Parallel Pin

K D Ø2.0	Unit mm Scale 1 : 1	
Code	PP 2022	PP 2023

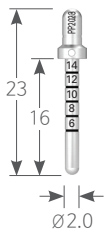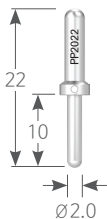

Surgical Tool

Drill

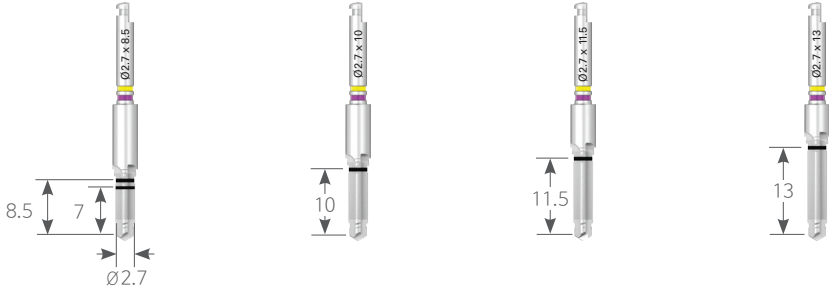K D Ø2.7		Unit mm Scale 1 : 1		
Code	DSU 2708M	DSU 2710M	DSU 2711M	DSU 2713M
				


Tap Drill

K		Unit mm Scale 1 : 1	
Fixture Size		Ø3.0	Ø3.3
Code		UNTD 3015	UNTD 3315
			

Fixture Driver

· Maximum Tightening Torque : 50Ncm

● Sold separately

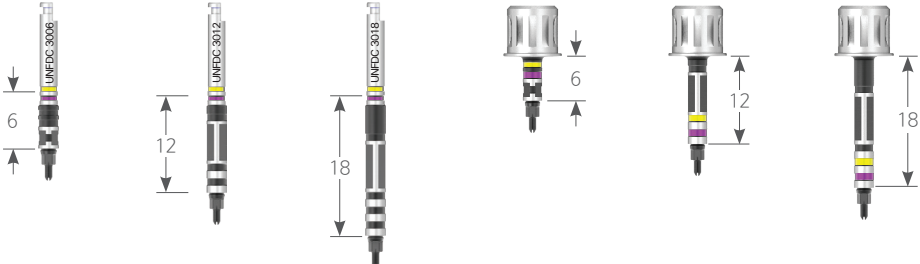K		Unit mm Scale 1 : 1			
Type		Machine		Ratchet	
Code	UNFDC 3006	UNFDC 3012	UNFDC 3018	UNFDW 3006	UNFDC 3012 UNFDC 3018
					

Screw Driver

K


Hex 1.2

Unit mm | Scale 1 : 1

Type	Machine	Ratchet
Code	MD 1230	HD 1220A



16.5

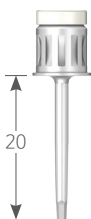

20

Solid Abutment Driver


K

D Ø4.0


Unit mm | Scale 1 : 1

● Sold separately

Length	6	12
Code	UNHDS 4006A	UNHDS 4012A

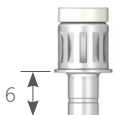

6

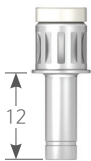

12

Surgical Protocol

Starter Drill

1,200rpm

Initial Drill

1,000rpm

Drill


800rpm

Tap Drill

50rpm

UF(II) Narrow Fixture

20rpm

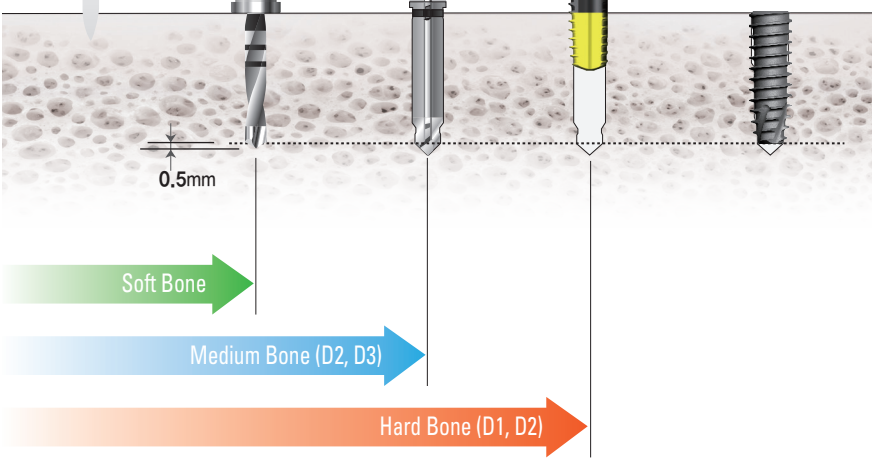

Tap Drill

1.5mm

4mm

· Tightening torque : 20~40Ncm


Maximum tightening torque : 50Ncm

UFII Wide Kit

Kit Code | UF 09

- Kit composed of drill with outstanding cutting force and durability, etc.
- UFII Wide Ø5.9 Ø6.4 Ø6.9 Fixture

Unit mm | Drill Scale 1 : 0.9

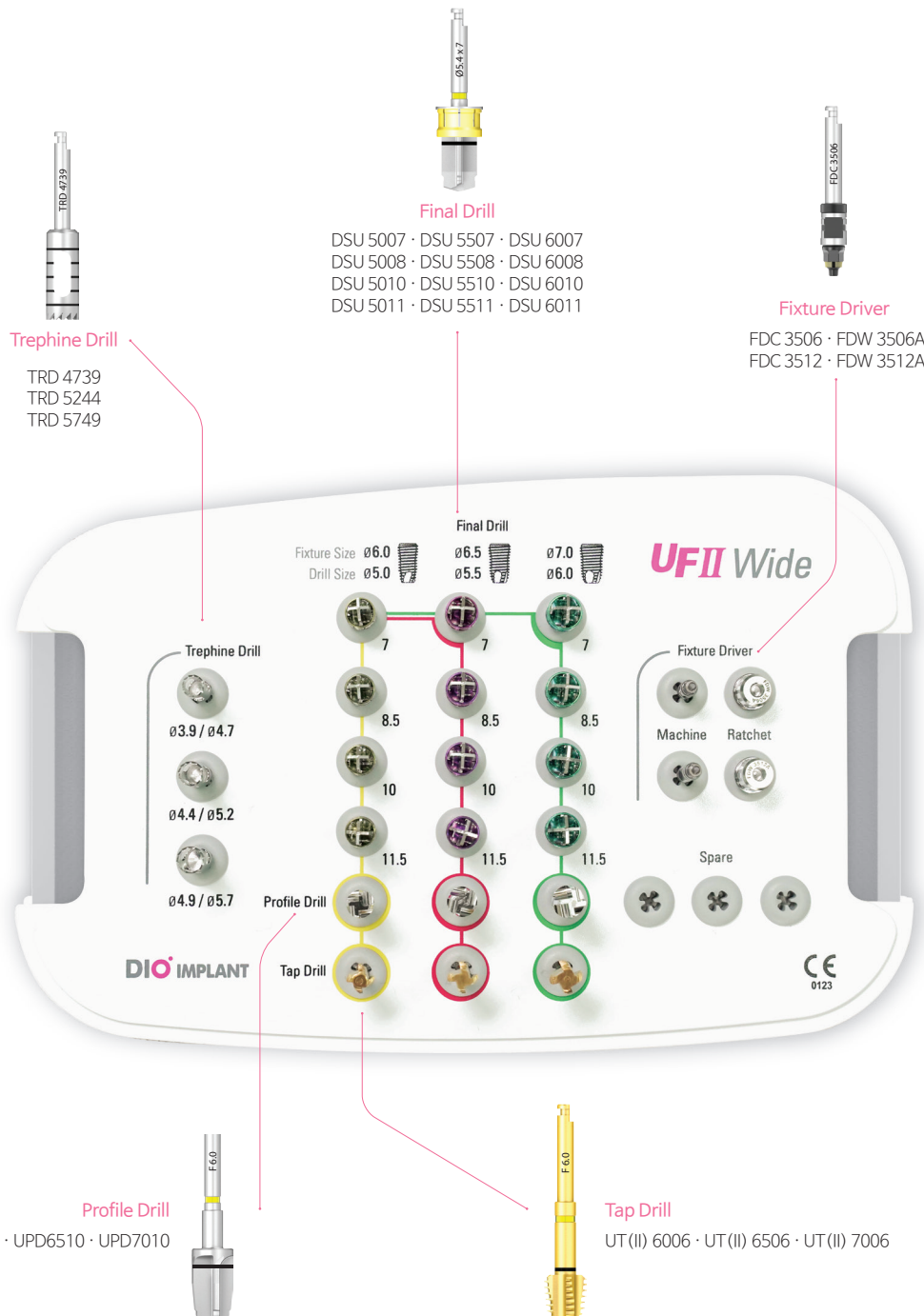

Surgical Tool

Trephine Drill



Scale 1 : 1.15

K		Unit mm Scale 1 : 1		
Code		TRD 4739	TRD 5244	TRD 5749
				
Diameter 1		Ø3.9	Ø4.4	Ø4.9
Diameter 2		Ø4.7	Ø5.2	Ø5.7

Profile Drill



Scale 1 : 1.15

K		Unit mm Scale 1 : 1		
Fixture Size		Ø6.0	Ø6.5	Ø7.0
Code		UPD 6010	UPD 6510	UPD 7010
				
Color		●	●	●

Surgical Tool

Final Drill



● Sold separately

K D Ø5.0						
Code	DSU 5007	DSU 5008	DSU 5010	DSU 5011	DSU 5013	DSU 5015
						
Length	7	8.5	10	11.5	13	15

K D Ø5.5						
Code	DSU 5507	DSU 5508	DSU 5510	DSU 5511	DSU 5513	DSU 5515
						
Length	7	8.5	10	11.5	13	15

K D Ø6.0						
Code	DSU 6007	DSU 6008	DSU 6010	DSU 6011	DSU 6013	DSU 6015
						
Length	7	8.5	10	11.5	13	15

Tap Drill

K

Unit mm | Scale 1 : 1

Fixture Size	Ø6.0	Ø6.5	Ø7.0
Code	UT(II) 6006	UT(II) 6506	UT(II) 7006

Color












Fixture Driver

K D Ø2.7

Unit mm | Scale 1 : 1

Type	Machine			Ratchet		
Code	FDC 3506	FDC 3512	FDC 3518	FDW 3506A	FDW 3512A	FDW 3518A

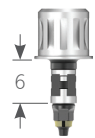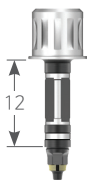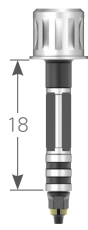

Surgical Protocol

Fixture Length 7mm

Extraction Socket &
Replace failed Implant

